


# 5-8-112 HOOKUP CONNECTOR ASSEMBLY REPLACEMENT INSTRUCTIONS

Using MEDIVATORS® ADVANTAGE and MEDIVATORS® ISA® Automated Endoscope Reprocessors

## HOOKUP DIAGRAM

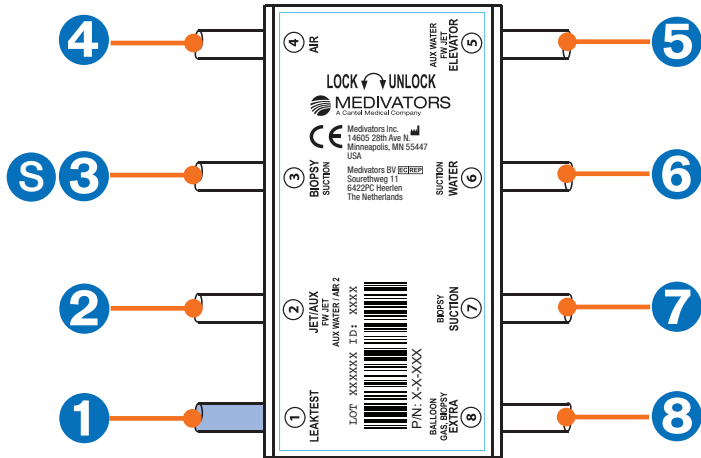

## CONNECTOR ASSEMBLY IDENTIFICATION CHART

Port ID	Port Name	Order Number	Connectors
1	LEAK TEST	H58112-1	
2	JET/AUX	N/A	N/A
3	BIOPSY	H58112-3	
4	AIR	N/A	N/A
5	ELEVATOR	N/A	N/A
6	WATER	N/A	N/A
7	BIOPSY	H58112-7	
8	EXTRA	N/A	N/A
S	N/A	N/A	N/A

If a Connector Assembly needs replacement, locate the port number/name and use the corresponding order number found in the chart above. Connector Assemblies are specific to individual hookups.

## DIRECTIONS FOR REPLACEMENT OF CONNECTOR ASSEMBLIES:

1. Identify if the Hookup Connector Block is the new (Figure 1 below) or previous (Figure 2 below) design.
2. Identify current Connector Assembly on the Hookup Connector Block to be replaced.
3. Remove current Connector Assembly by pulling associated tubing **B** from the Hookup Connector Block Barb **A** (NOTE: If using the previous design, unscrew the Connector Coupling **F** of Connector Assembly from the Port Connector **C**, keeping the Port Connector firmly held in place and pull the tubing off the Port Connector stem **D**). Discard this Connector Assembly.
4. Retrieve the enclosed new Connector Assembly. Note that it may not be identical to the Connector Assembly it is replacing. Discard Ferrule **E** and Connector Coupling **F** if using the new design. Confirm Connector Assembly labeling matches the Hookup Connector Block Label Port Name and the Port ID Number using the Connector Assembly Identification Chart. NOTE: Port ID # may not be identified on some older Hookup Connector Block labels.
5. Push tubing **B** of the new Connector Assembly onto the Hookup Connector Block Barb **A**. NOTE: If experiencing difficulty in attaching the tubing to the Hookup Connector Block Barb **A**, dip the tip of the tubing into hot water for 10 seconds. (NOTE: If using the previous design, push tubing of new Connector Assembly onto the Port Connector Stem **C** and manually screw Connector Coupling **F** onto Port Connector **C** until tight. Do not over tighten. Ensure Port Connector **C** does not turn. If needed, use a wrench to prevent turning.)

## HOOKUP CONNECTOR BLOCK TYPES

### Exploded View of Port Connector, Coupling Connector, Ferrule and Tubing:

Figure 1 - New Design


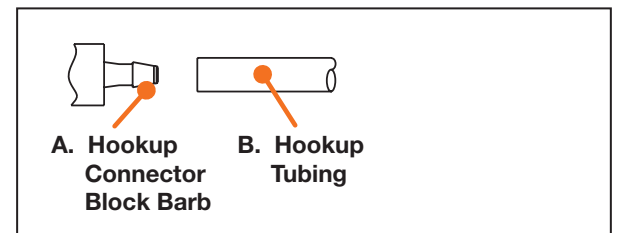


Figure 2 - Previous Design



## WARNINGS & CAUTIONS

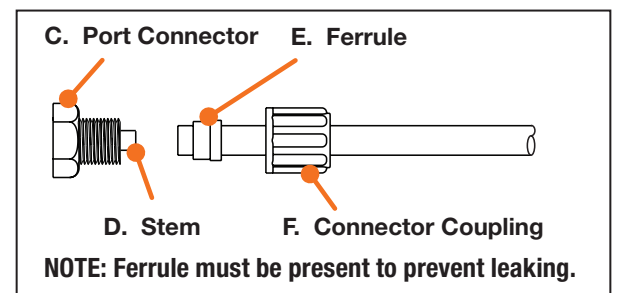

This hookup is intended for connection only with the specific endoscope models identified in the appropriate MEDIVATORS online hookup guide. Use of this hookup with endoscopes other than those specified by the online hookup guide may result in inadequate HLD.

Users operating MEDIVATORS AERs and hookups must be trained and competent in the understanding of endoscope channel systems.

Prior to reprocessing, users must verify that the hookup used contains connections for all appropriate channels that require a separate reprocessing connection according to the hookup instructions.

All connectors and adapters must remain firmly attached and unrestricted for the entire disinfection cycle to ensure adequate HLD. Users must inspect the endoscope and hookup combination to verify proper connection and flow.

Modifications or repairs to hookups which do not correspond to the manufacturer's specifications may result in inadequate HLD and/or damage to endoscopes.

Failure to properly clean and prepare endoscope for immersion prior to hookup connection may result in inadequate HLD and/or damage to the endoscope.

Endoscopes must be inspected for damage and verified to be in proper working order prior to hookup connection. Connection of hookups to damaged endoscopes or endoscopes with obstructed/restricted channels could result in inadequate HLD.

Manual leak testing must be carried out in accordance with manufacturer's recommendations prior to reprocessing.

Install all necessary waterproof caps, plugs and cleaning adapters prior to immersion.

Do NOT autoclave MEDIVATORS hookups.

MEDIVATORS makes no claim on high-level disinfection efficacy when these reprocessing instructions are not followed, or when this hookup kit is applied to endoscopes other than those specified in the online interactive hookup guide.

Avoid biological contamination. Always wear personal protective equipment when handling endoscopes.

If you have any questions regarding this hookup, in the U.S., please contact Medivators Technical Support at 1-800-444-4729. Outside the U.S., please contact your Medivators representative.

MEDIVATORS® is a registered trademark of Medivators Inc.  
ISA® is a registered trademark of Cantel Medical (Italy) S.r.l.  
All company or product names referenced are trademarks of their respective owners.



Medivators Inc.  
14605 28th Avenue North  
Minneapolis, MN  
55447-4822 USA  
Tel: +1.763.553.3300  
Fax: +1.763.553.3387

Medivators BV  
Sourethweg 11  
6422 PC Heerlen  
The Netherlands  
Tel: +31.45.5.471.471  
Fax: +31.45.5.429.595

Cantel (France) SAS  
3 rue du pré fermé  
31200 Toulouse  
France  
Tel: +33 (0) 5.61.245.275  
Fax: +33 (0) 5.61.247.001

Cantel Medical (Italy) S.r.l.  
Via Laurentina, 169  
00071 Pomezia (RM) Italia  
Tel.: +39 06 9145399  
Fax: +39 06 9146099  
www.cantelmedical.it

Cantel Medical (UK) Limited  
Wolsley House  
Staffordshire Technology Park  
Beaconside, Stafford ST18 0GA  
United Kingdom  
Tel: +44 (0) 1785.782420  
Fax: +44 (0) 1785.782421

Cantel Medical Asia/Pacific Pte. Ltd.  
1A International Business Park  
#05-01 Singapore 609933  
Tel: +65.6227.9698  
Fax: +65.6225.6848

Cantel Medical Devices (China) Co. Ltd.  
Unit 804-805, Innov Tower Block A,  
Hongmei Road, Xuhui 200233 Shanghai  
Tel: +86 21 60161380  
Fax: +86 21 61210913

www.medivators.eu