

Automated, safe disinfectant neutralization



Endoscope reprocessing the way it should be

MEDIVATORS



Provides safe containment, neutralization

Environmentally responsible disinfectant neutralization system

- Safe - protects staff by containing disinfectant fumes
- Flexible – neutralizes OPA and Glutaraldehyde disinfectants
- Fast – neutralizes in as little as 5 minutes
- Easy to use, semi-automated operation
- Compliance - meets state and local neutralization regulations

Neutra-Hyde™ chemical neutralizer



Ready-to-use, easy-pour neutralizer bottle



ion and disposal of used disinfectants

Step-by-step instructions and large control switches for ease of operation

Oversized cap for easy adding of neutralizer


Large capacity 12-gallon neutralization tank

Built-in shelf holds Neutra-Hyde neutralizer

Quick connect fittings provide safe and secure disinfectant transfer

Sealed reservoir contains disinfect fumes

Compact design for ease of mobility

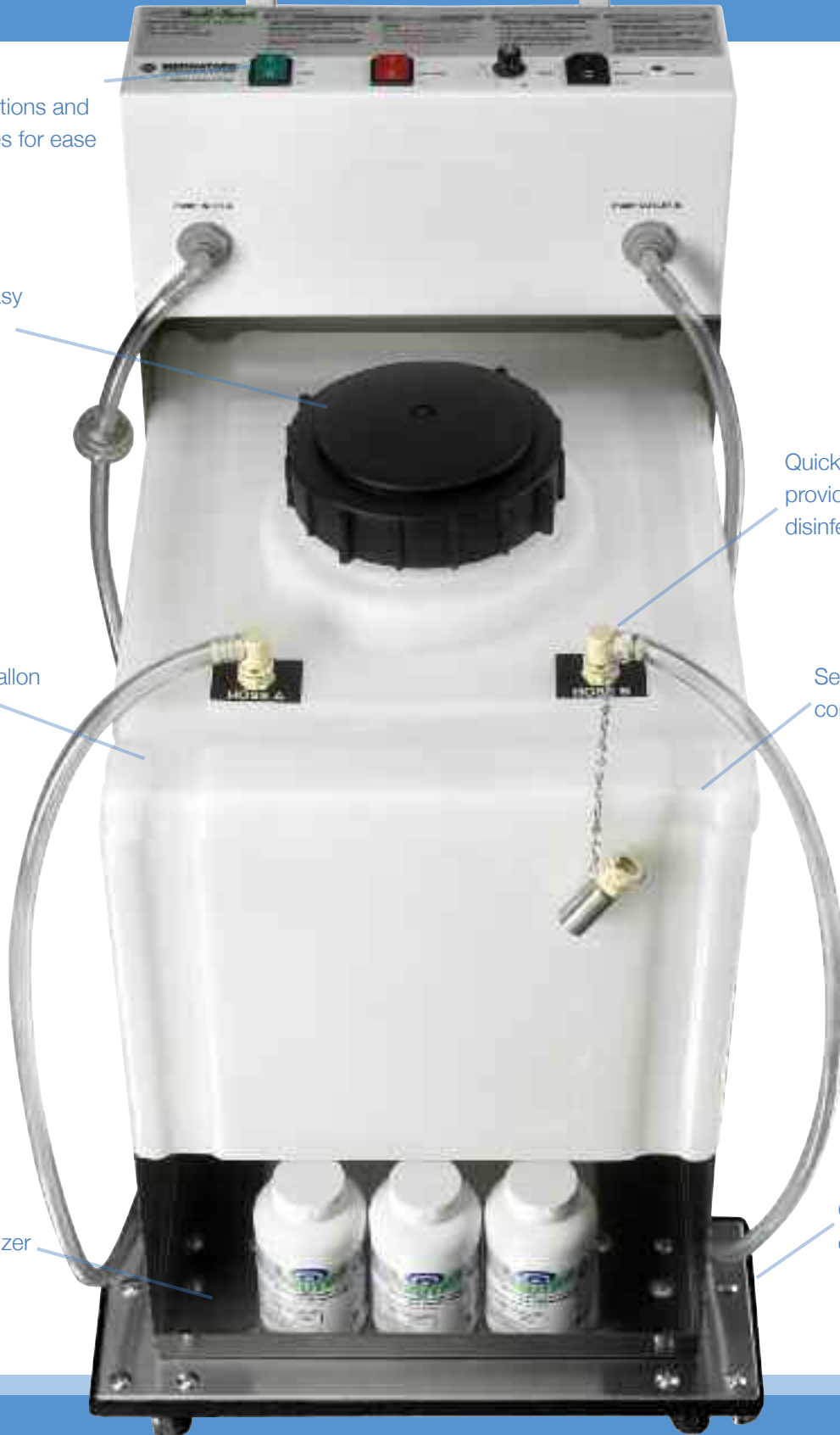



17.5" (44.5 cm)

Dimensions

37 H x 17.5 W x 20.75 D (inches)
94 H x 44.5 W x 52.7 D (cm)

Weight

35 lbs (15.8 kg)

Electrical Requirements

100-120VAC, 50/60Hz, 0.8 amp

Electrical Safety Certifications

UL Standard 61010-1
CSA C22.2 Standard 61010-1
NRTL Listed
MET Labs listing #E112640

Unit Placement

Mobile cart

Tank Capacity

12 gallons (45.3 liters)

Warranty

One year



20.75" (52.7 cm)

37"
(94 cm)



Medi-Newt Ordering Information

Part number	Description
DNS-1001	Medi-Newt System
DNS-1002	Neutra-Hyde Neutralizer (12 bottles per case)
DNS-1003	Tubing Kit (hose A and hose B)
DNS-1004	CER Reservoir Adapter Tubing
DNS-1005	In-Line Filters (5 per pack)



Minntech Corporation
14605 28th Avenue North
Minneapolis MN 55447
Tel: +1-763-553-3300
Fax: +1-763-553-3387
www.minntech.com

Minntech Asia/Pacific Pte Ltd
438 Alexandra Road
#03-01 Alexandra Point
Singapore 119958
Tel: +65-6227-9698
Fax: +65-6225-6848

Minntech BV
Sourethweg 11
6422 PC Heerlen
The Netherlands
Tel: +31-45-5-471-471
Fax: +31-45-5-429-695

Minntech China, Rep. Office
Room 1508B, Avic Plaza
No. 10 Yi Dongsanhuan Zhonglu
Chaoyang District,
Beijing 100022 China
Tel: +86-10-6567-8446
Fax: +86-10-6567-8445

For additional information please contact your Medivators representative or visit our website at www.minntech.com/medivators