

# ADVANTAGE<sup>®</sup> PLUS



45" (114 cm)  
49" (124 cm) w/barcode bracket

## Electrical Supply

120 VAC ± 10%: single phase, 60HZ; 15A (North America)  
240 VAC ± 10%: single phase, 50/60HZ; 10A (International)

## Water Supply

Minimum flow rate: 1.32 GPM (5 LPM) with a pressure of 35 to 87 PSI (2.4 to 6 bar)  
Water Temperature: 95°±4°F (35° ± 2°C)  
Hardness: less than 12 gpg (200 ppm)  
Pre-filter wall space: 26 H x 38 W x 6 D inch (66 H x 96 W x 15 D cm) within 3 feet (100 cm) of the Advantage Plus

## Drain

3.17GPM (12 LPM) less than 18" (46 cm) from the floor, vented

## Air Supply

Compressed air (dry/oil free): 58 to 125 PSI (4 to 8.6 bar), minimum flow volume of 1.5 CFM (0.88 m<sup>3</sup>/h)  
Optional air compressor, adds 11" (28 cm) to depth of unit

## Basin Size


19 W x 17 L x 7 D inch (48 H x 43 L x 18 D cm)

## Advantage Plus Ordering Information - North America

| Part number | Description   |
|-------------|---|
| ADV-1007    | Advantage Plus for Rapicide <sup>®</sup> PA with air compressor |
| ADV-1008    | Advantage Plus for Rapicide PA without air compressor           |
| ML02-0116   | Rapicide PA Disinfectant (Canada)                               |
| ML02-0117   | Rapicide PA Disinfectant  |
| ML02-0118   | Rapicide <sup>®</sup> PA test strips                            |

## Advantage Plus Ordering Information - International

| Part number | Description   |
|-------------|---|
| ADV-2001    | Advantage Plus for Rapicide PA with air compressor          |
| ADV-2002    | Advantage Plus for Rapicide PA without air compressor       |
| ADV-4001    | Advantage Plus for Adaspor <sup>®</sup> with air compressor |
| ADV-4002    | Advantage Plus for Adaspor without air compressor           |
| ML02-0115   | Rapicide PA Disinfectant (International)                    |
| ML02-0118   | Rapicide PA test strips                                     |

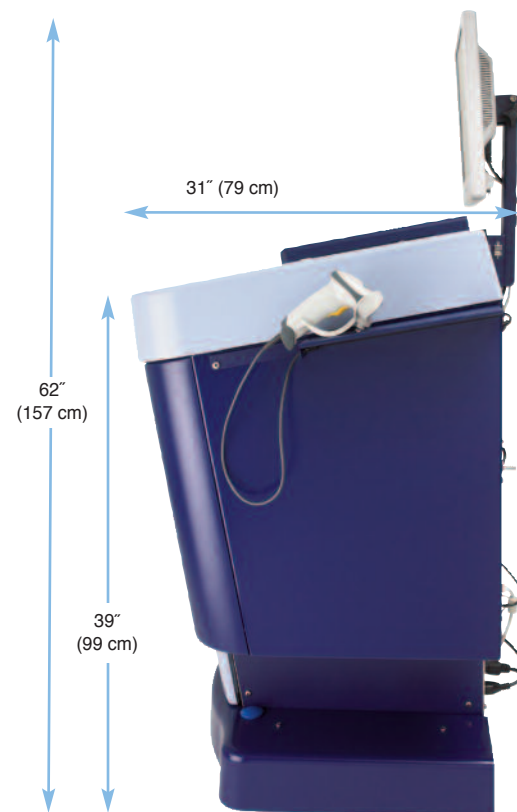

# ADVANTAGE<sup>®</sup> PLUS



**MINNTECH<sup>®</sup>**  
A Cantel Medical Company

Minntech Corporation  
14605 28th Avenue North  
Minneapolis MN 55447  
Tel: +1-763-553-3300  
Fax: +1-763-553-3387

Minntech Asia/Pacific Pte Ltd  
No.1 A International Business Park  
#05-01 Singapore 609933  
Tel: +65-6227-9698  
Fax: +65-6225-6848

Minntech BV  
Sourethweg 11  
6422 PC Heerlen  
The Netherlands  
Tel: +31-45-5-471-471  
Fax: +31-45-5-429-695

Minntech China, Rep. Office  
Room 708, 7th Floor Kaiheng Center,  
No. 2 Chaoyangmennei Street,  
Dongcheng District, Beijing China, 100010  
Tel: +8610-6567-8446  
Fax: +8610-6567-8445

For additional information please contact your Medivators representative  
or visit our website at [www.minntech.com/medivators](http://www.minntech.com/medivators)

Endoscope reprocessing the way it should be

**MEDIVATORS<sup>®</sup>**

# ADVANTAGE<sup>®</sup> PLUS

The most advanced Automated Endoscope Reprocessor from the market leader

## Rapicide<sup>®</sup> PA High-Level Disinfectant

- Proven endoscope compatibility
- Faster endoscope availability
- Aldehyde-free
- 30°C disinfection
- Reduced worker exposure
- Environmentally friendly



## Large, easy loading, independent basins

- Asynchronous operation
- Accommodates all flexible endoscopes

## Automated alcohol flush



- Facilitates drying of endoscope channels
- Reduces incidence of microbial growth

## Automated leak testing

- Enables disinfection of an endoscope prior to repair
- Minimizes leak-related repair costs

## High performance cleaning cycle

- Enhances manual cleaning and optimizes disinfection

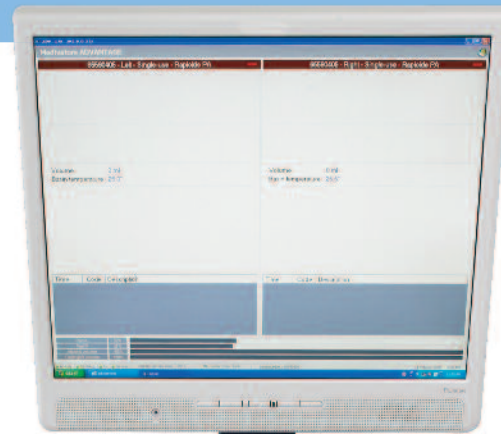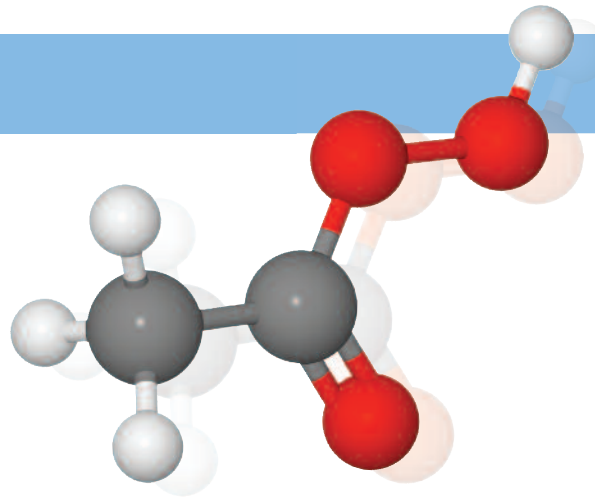

## Comprehensive endoscope channel monitoring


- Guarantees proper hookup connection
- Monitors endoscope channels for blockage
- Ensures all channels are correctly and fully reprocessed

## Extensive cycle data management

- Barcode reader for fast, accurate data entry
- PC-based for scope-specific cycle programming, operation, and cycle data documentation
- Real-time, on-screen cycle monitoring

## Networking & remote diagnostics

- Network multiple systems together
- Connect to local or off-site server for centralized data collection
- Fast technical support via on-line diagnostics

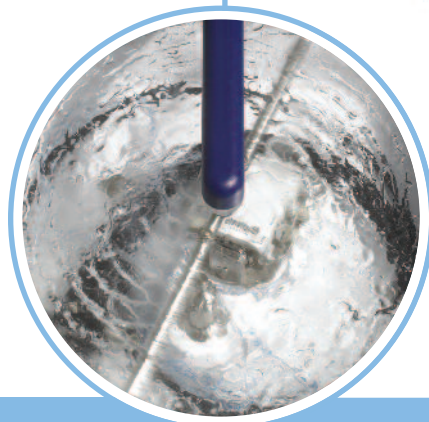

## Hands-free lid operation

- Provides ease of use
- Lessens risk of cross contamination